Boppin' Beaver Switch Adapted Toy USER QUICK GUIDE




Introduction

The switch adapted Boppin' Beaver toy sings and dances and can be controlled with accessible switches that have a 3.5 mm jack. The beaver is a Fisher Price toy that teaches ABCs, 123, opposites, and simple games.

Features

Power button

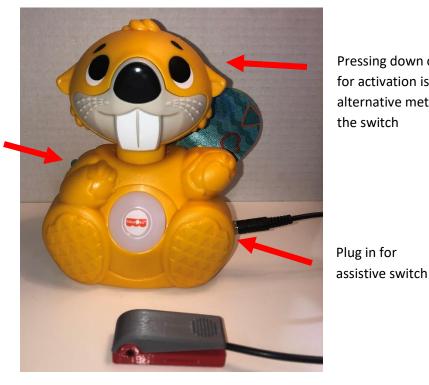

Pressing down on the head for activation is the alternative method to using the switch

Plug in for

Usage

To activate the beaver, the green power switch must be turned on, which is located on the beaver's shoulder. Once turned on, the assistive switch can be used to get the beaver to do the ABCs, 123, opposites, and simple games.

Compatibility

The Boppin' Beaver switch adapted toy is compatible with any switches that have 3.5 mm cable ends. The beaver can be powered with three AA batteries.

V1.0 | APRIL 2022

Boppin' Beaver Switch Adapted Toy USER QUICK GUIDE



Specifications

Height (mm)	Width (mm)	Depth (mm)	Weight (g)
145	110	80	413

Cleaning

Wipe the beaver with a clean cloth. Be careful around the exposed electronic plug ins if the cloth is damp.